




CHEMGRATOR


Integrator for Steam Sterilization

Ref.: 134 ° C

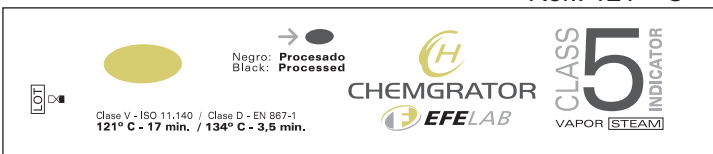

0 Minutes

Ref.: 121 ° C

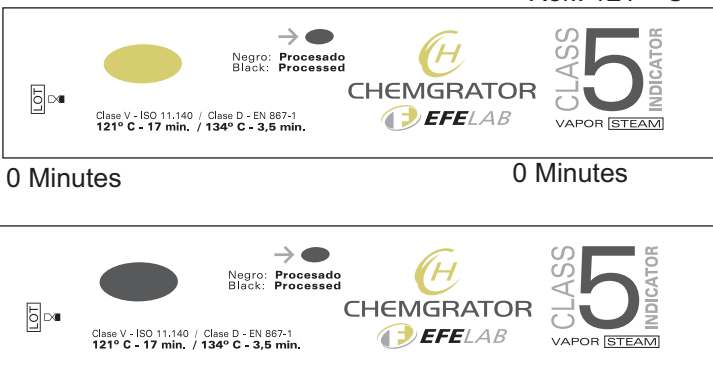

0 Minutes

3.5 Minutes

17.5 Minutes

CHEMGRATOR integrating indicators are specially designed to integrate the effects of time, temperature and saturated steam in sterilization processes.

The change of color of the integrator occurs in the presence of these three conditions providing visual evidence immediately after if the sterilization conditions have been achieved during the cycle. This can be easily seen by the change of color of the indicator up to its reference color

Main features of CHEMGRATOR steam integrating indicators:

Hi-quality indicator ink that meets the requirements for ISO 11140-1, EN-867 and IRAM 37101-1 Class 5 standards.

An especially formulated ink changes from green to black and the ending color is printed on each integrator for an immediate interpretation of the results of the cycle.

The color change is irreversible. Thus it can be used as a record of the sterilization cycle.

Suitable for regular and flash steam sterilization cycles.

To be used in a steam autoclave at 121° C - 17 minutes or at 134° C - 3.5 minutes

Shelf Life: 60 months from its date of manufacture.